

BELOURA

BODY FILLER



BELOURA
Premium
Body Contouring Filler

Type	Cross-linked Hyaluronic Acid Filler	Effect	Enhances and restores natural body volume
Hyaluronic Acid Concentration	20mg/ml	Indicated Areas	Indicated for breast, buttocks, and soft-tissue contouring
Volume	60 ml / vial	Shelf Life	36 months
Appearance	Transparent monophasic gel	Needle Size	16~18G
Lidocaine Content	3mg/ml		

**Recommended
For**



Those who desire fuller,
more defined body contours



Those who wish to improve
breast firmness and elasticity



Those seeking quick and
natural volume enhancement



For those who prefer a safe,
non-surgical solution

Premium Body Filler for Natural Elegance
BELOURA — Key Advantages



Safety



Longer Duration



Higher Viscosity



High Quality



Body Contouring
Defines
natural body lines

Hips Sculpting
Enhances smooth
and lifted curves

Skin Hydration
Improves texture
and elasticity


Breast Volumizing
Restores
natural fullness

Buttocks Augmentation
Creates attractive
body balance

Why BELOURA?

BELOURA is manufactured under the strict supervision of the **Korean MFDS** and certified with **ISO 9001** and **ISO 13485**.

BELOURA Manufacturing Process

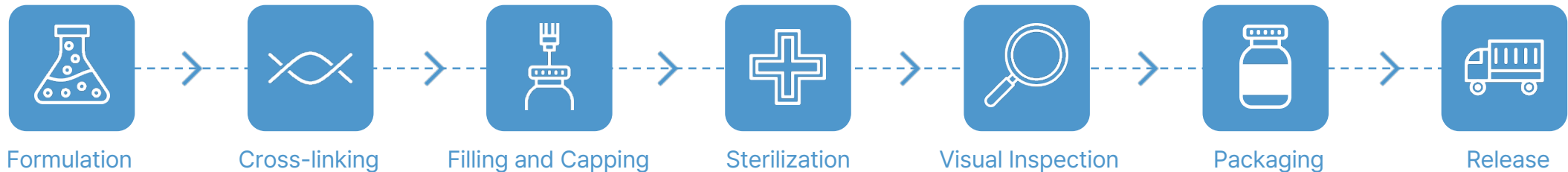

Good Quality Raw Material

BELOURA contains **high-purity, non-animal Hyaluronic Acid**, registered with **EDQM**, and formulated using **WFI (Water for Injection)** instead of conventional purified water — ensuring exceptional safety and purity.

BELOURA features a **monophasic gel structure**, delivering **exceptional results** for deep wrinkle correction and body contouring.

With a **storage modulus (G')** that is **over three times higher** than ordinary fillers, it **degrades gradually** and provides **long-lasting volume, elasticity, and natural beauty**.

BELOURA is manufactured using **innovative DLT and UDL technologies** with an advanced **purification process** that significantly reduces impurities and particle size. This enhances **purity and injectability**, while minimizing the risk of side effects and ensuring **stable product quality**.

DLT: Dialysis at Low Temperature

UDL: Ultra-fine Dispersion at Low Temperature

Simple Preparation Steps



1 Prepare one vial of BELOURA, a 10 ml syringe, saline solution, and a 16–18 G cannula.



2 Draw 10 ml of sterile saline into the syringe.



3 Inject the saline into the BELOURA vial and mix gently.




4 The filler is now ready for smooth extraction and use.

Usage Precautions

1. Ensure the treatment area is thoroughly disinfected before the procedure.
2. The product may be diluted with sterile saline for smoother extraction, but the entire process must be performed under strict aseptic conditions.
3. Inject only into the subcutaneous fat layer of the desired area, such as the hips. Do not inject into blood vessels under any circumstances.
4. When injecting into the breast, avoid any contact with the mammary glands.
5. Recommended dosage per area:
Minimum 50cc – Maximum 200cc for each breast or hip.
6. After the procedure, maintain strict hygiene and infection prevention care to ensure optimal healing.

Injection Skin Layer

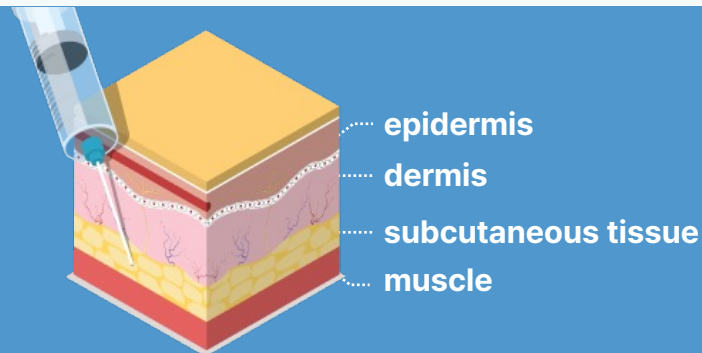

BELOURA is designed for injection into the subcutaneous layer to enhance natural body volume.

Before & After

Breast Volumizing



Before

After

Buttocks Contouring



Before

After



KFLUX Co.,Ltd

Room 909, Woorim e-Biz Center 1, 28 Digital-ro 33-gil, Guro-gu, Seoul,
Republic of Korea

☎ +82 10-4919-5228 ✉ kflux@naver.com 🏠 www.kflux.co.kr

WhatsApp:

